# **Human CTDSP2 Knockdown Cell Line (WB-Validated)**



**Catalog #: C63822** 

#### **Aliases**

CTDSP2; CTD Small Phosphatase 2; SCP2; OS4; Nuclear LIM Interactor-Interacting Factor 2; NLI-Interacting Factor 2; Small CTD Phosphatase 2; PSR2; CTD (Carboxy-Terminal Domain, RNA Polymerase II, Polypeptide A) Small Phosphatase 2; Carboxy-Terminal Domain RNA Polymerase II Polypeptide A Small Phosphatase 2; Conserved Gene Amplified In Osteosarcoma; Small C-Terminal Domain Phosphatase 2; Protein OS-4; EC 3.1.3.16; NIF2

## **Background**

Gene Name: CTDSP2 NCBI Gene Entry: 10106

## **Storage**

Store at liquid nitrogen for 1 year.

## **Kit Components**

- 1. Human CTDSP2 Knockdown Cell Line (Wb-Validated)
- 2. Wild-type cell line

### **Parental Cell Line**

Human cell line supplied by the client

#### **Validation Methods**

RT-qPCR, Western blotting (WB)

## **Shipping**

Shipped on Dry Ice.


### **Instructions For Use**

This knockdown cell line should be paired with wild-type cell line for use.

**Note:** This product is for research use only.


#### Validation Data

## **Human CTDSP2 Knockdown Cell Line (WB-Validated)**

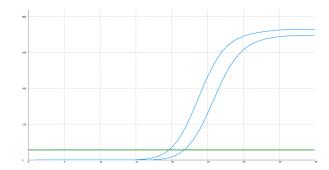

Genotype	Ct Value
Wild-Type	19.47
Knock-Down	21.50
$\Delta Ct (Ct_{KD}-Ct_{WT})$	2.03
% mRNA	
Reduction	<b>.</b> 76%

RT-qPCR analysis. HeLa cells were infected with CTDSP2-specific shRNA lentiviral particles, total RNA was extracted from wild-type and knockdown cells, RT-qPCR was performed using gene-specific primers.  $\Delta$ Ct (CtKD-CtWT) was used to calculate mRNA reduction (%) between wild-type and knockdown cells using the following formula:  $(1-1/2\Delta$ Ct) x 100%.

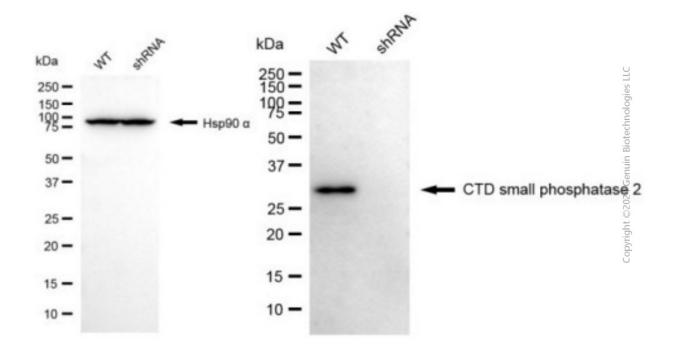

Western blotting analysis. CTDSP2 protein expression in wild-type (WT) and shRNA knockdown (KD) HeLa cells was detected using Western blotting. Hsp90  $\alpha$  served as a loading control. The blots were incubated with primary antibodies against CTDSP2 and Hsp90  $\alpha$ , respectively, followed by incubating with HRP-conjugated goat anti-mouse secondary antibody. Images were developed using FeQ<sup>TM</sup> ECL Substrate Kit.