Human NFKBIE Knockdown Cell Line (WB-Validated)



Catalog #: C65670

Aliases

NFKBIE; NFKB Inhibitor Epsilon; IKBE; NF-Kappa-B Inhibitor Epsilon; Nuclear Factor Of Kappa Light Polypeptide Gene Enhancer In B-Cells Inhibitor, Epsilon; I-Kappa-B-Epsilon; IkappaBepsilon; NF-Kappa-BIE; IkB-Epsilon; IkB-E

Background

Gene Name: NFKBIE NCBI Gene Entry: 4794

Storage

Store at liquid nitrogen for 1 year.

Kit Components

- 1. Human NFKBIE Knockdown Cell Line (Wb-Validated)
- 2. Wild-type cell line

Parental Cell Line

Human cell line supplied by the client

Validation Methods

RT-qPCR, Western blotting (WB)

Shipping

Shipped on Dry Ice.


Instructions For Use

This knockdown cell line should be paired with wild-type cell line for use.


Note: This product is for research use only.

Validation Data

Human NFKBIE Knockdown Cell Line (WB-Validated)

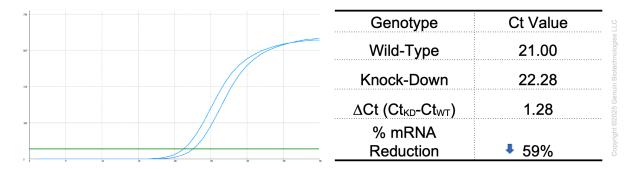

RT-qPCR analysis. HeLa cells were infected with NFKBIE-specific shRNA lentiviral particles, total RNA was extracted from wild-type and knockdown cells, RT-qPCR was performed using gene-specific primers. Δ Ct (CtKD-CtWT) was used to calculate mRNA reduction (%) between wild-type and knockdown cells using the following formula: $(1-1/2\Delta$ Ct) x 100%.

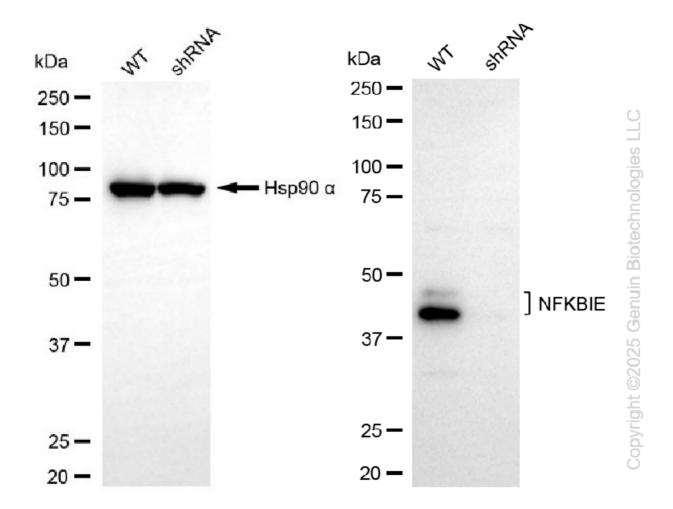

Western blotting analysis.NFKBIE protein expression in wild-type (WT) and shRNA knockdown (KD) HeLa cells was detected using Western blotting. Hsp90 α served as a loading control. The blots were incubated with primary antibodies against NFKBIE and Hsp90 α , respectively, followed by incubating with HRP-conjugated goat anti-rabbit secondary antibody. Images were

TEL: +1-540-855-7041

PAGE 3

Human NFKBIE Knockdown Cell Line (WB-Validated)

developed using FeQ $^{\text{TM}}$ ECL Substrate Kit.