WB-Validated ISG15 Knockdown Cell Lysate Kit



Catalog #: L62197

Aliases

ISG15; ISG15 Ubiquitin Like Modifier; UCRP; IFI15; G1P2; Interferon, Alpha-Inducible Protein (Clone IFI-15K); Ubiquitin Cross-Reactive Protein; Ubiquitin-Like ProteinISG15; HUCRP; IP17; Interferon-Induced 17-KDa/15-KDa Protein; Interferon-Stimulated Protein, 15 KDa; Interferon-Induced 15 KDa Protein; Interferon-Induced 17 KDa Protein; IMD38

Background

Gene Name: ISG15 NCBI Gene Entry: 9636

Storage

Stored at -20°C for 2 years.

Kit Components

1. 100 µg WT cell lysate

2. 100 µg KD cell lysate

Tested Cell Line

HepG2

Validation Methods

RT-qPCR; Western Blotting (WB)

Shipping

Shipped with gel ice packs. Immediately store the product in a standard freezer at -20°C upon receipt.

Instructions For Use

This knockdown cell lysate should be paired with wild-type HepG2 cell lysate for use. For Western blotting, we recommend running wild-type and knockdown lysates on the same gel, and loading each well with equal volume and equal amount of total proteins.

Manufacturing Process

The following protocol was used to generate mRNA knockdown cell lysate:

1.Release 0.5 million HepG2 cells into a 35 mm tissue culture dish in 2 mL of the growth medium


WB-Validated ISG15 Knockdown Cell Lysate Kit

- (DMEM containing 10% FBS and 1% pen/strep). Cell density should reach 50-60% confluence the following day.
- 2.24 h after cell release, pre-warm the shRNA lentiviral medium to 37°C.
- 3.Discard 1 mL of the original growth medium of the 35 mm dish.
- 4. Using a serological pipette, gently mix the lentiviral solution 3 times.
- 5.Carefully add 1 mL of the lentiviral solution to the well. Tip: To prevent splashing, add the solution to the dish along the wall.
- 6.Add a polybrene stock solution to the culture medium at a final concentration of 5 $\mu g/mL$. Gently swirl the dish to mix.
- 7.48 h after cell release, without discarding the original medium, add another 1 mL of lentiviral medium directly into the dish.
- 8.Add an additional polybrene stock solution into the dish to obtain a final concentration of 5 μ g/mL. Tip: Now, the medium in the dish should be a total of 3 mL.
- 9.72 h after cell release, cells may reach confluence. Trypsinize the cells off the 35 mm dish and culture those cells in a 60 mm dish.
- 10.Add puromycin to the dish at a final concentration of 4 μ g/mL. Tip: To assess the efficacy of puromycin selection, culture a dish of wild-type HepG2 cells as a negative control.
- 11.Allow puromycin selection for 48 h. Almost all wild-type HepG2 cells should die, while the dish infected with lentiviruses should have some remaining cells.
- 12.Replace the medium with regular growth medium without puromycin and allow the cells to grow to confluence before harvesting or staining.

Cells were lysed with IntactProteinTM cell/tissue lysis kit (Cat#415) and stored in -20°C.

Note: This product is for research use only.


Validation Data

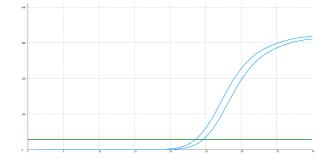

Genotype	Ct Value
Wild-Type	23.20
Knock-Down	24.23
∆Ct (Ct _{KD} -Ct _{WT})	1.03
% mRNA	
Reduction	J 51%

RT-qPCR analysis. HepG2 cells were infected with ISG15-specific shRNA lentiviral particles, total RNA was extracted from wild-type and knockdown cells, RT-qPCR was performed using gene-specific primers. Δ Ct (CtKD-CtWT) was used to calculate mRNA reduction (%) between wild-type and knockdown cells using the following formula: $(1-1/2\Delta$ Ct) x 100%.

WB-Validated ISG15 Knockdown Cell Lysate Kit

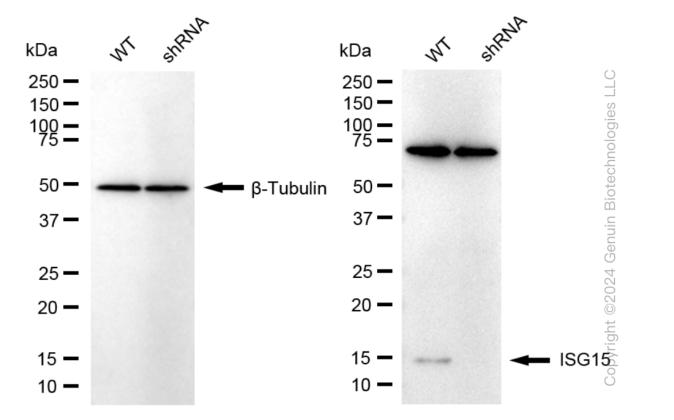

Western blotting analysis. ISG15 protein expression in wild-type (WT) and shRNA knockdown (KD) HepG2 cells was detected using Western blotting. β -Tubulin served as a loading control. The blots were incubated with primary antibodies against ISG15 and β -Tubulin, respectively, followed by incubating with HRP-conjugated goat anti-rabbit secondary antibody. Images were developed using FeQTM ECL Substrate Kit.