

## Catalog #: T013

### Aliases

HA-tag; HA epitope tag; HA1; HA2; hemagglutinin; Hemagglutinin HA1 chain; Hemagglutinin HA2

### Background

The HA-tag is a small, widely used epitope tag derived from the hemagglutinin (HA) glycoprotein of the human influenza virus. The most common HA-tag sequence (YPYDVPDYA, 9 amino acids) corresponds to an immunodominant epitope recognized by monoclonal antibodies. It is typically fused to the N-terminus or C-terminus of recombinant proteins via genetic engineering. Unlike larger fusion tags (e.g., GFP), its compact size minimizes interference with protein structure or function while enabling highly specific detection.

### Application Information

Molecular Weight: Recombinant protein dependent

Clonality: Mouse monoclonal antibody

Species Reactivity: Recombinant protein

Applications Tested: Western Blotting (WB), immunocytochemistry (IC), immunoprecipitation (IP)

### Immunogen

A synthesized peptide derived from HA-tag

### Storage Buffer

Supplied in PBS (pH 7.3) containing 30% glycerol, and 0.01% sodium azide.

### Storage

Store at -20 °C for one year.

### Recommended Dilutions

Western Blotting (WB): 1:2,000-1:5,000

Immunocytochemistry (IC): 1:200-1:500

Immunoprecipitation (IP): 1:100-1:200

**Note:** This product is for research use only.

### Validation Data

---


#### SUPPORT

SUPPORT@GENUINBIOTECH.COM  
TEL: +1-540-855-7041

#### ORDERS

SALES@GENUINBIOTECH.COM  
FAX: +1-540-855-7041

[WWW.GENUINBIOTECH.COM](http://WWW.GENUINBIOTECH.COM)

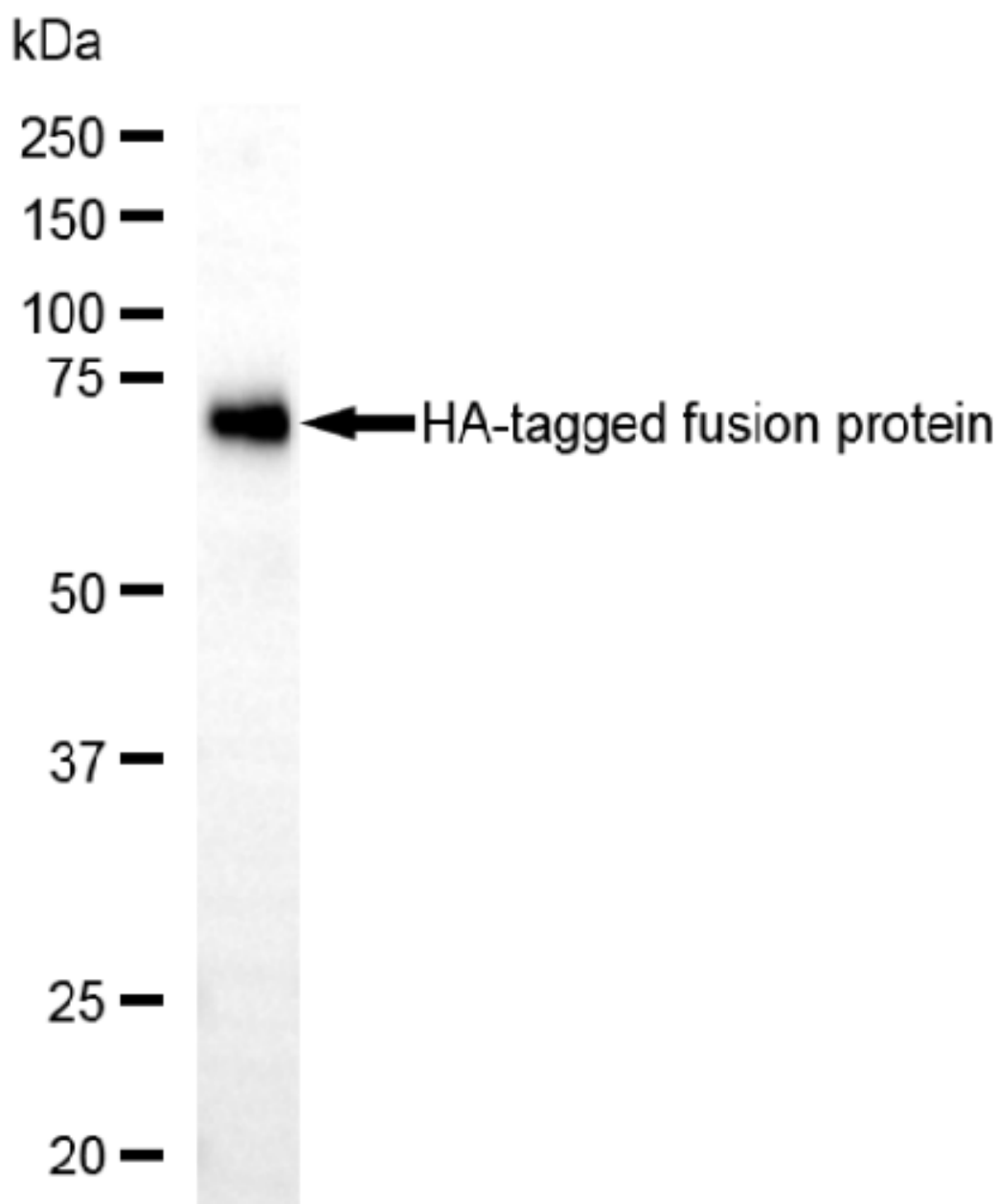

Copyright ©2025 Genuin Biotechnologies LLC

Western blotting analysis using anti-HA-tag antibody (Cat#T013). Recombinant multi-tag protein (5 ng) in E. coli were loaded and separated by SDS-PAGE. The blot was incubated with anti-HA-tag antibody (Cat#T013, 1:2,500) and HRP-conjugated goat anti-rabbit secondary antibody (Cat#201, 1:20,000) respectively. Image was developed using FeQ™ ECL Substrate Kit (Cat#226).


**SUPPORT**

SUPPORT@GENUINBIOTECH.COM  
TEL: +1-540-855-7041

**ORDERS**

SALES@GENUINBIOTECH.COM  
FAX: +1-540-855-7041

**WWW.GENUINBIOTECH.COM**

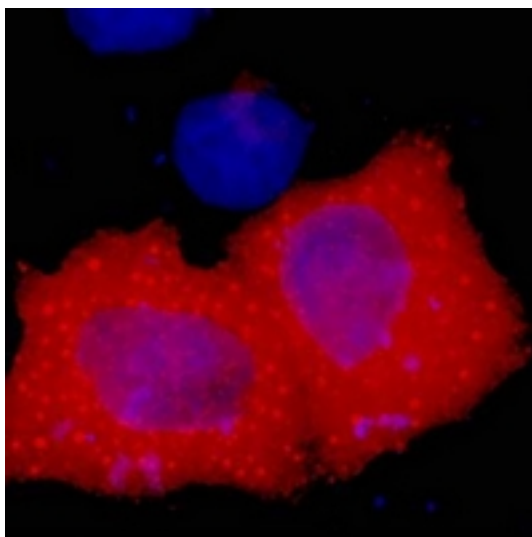

Immunocytochemical analysis of HA-tag staining in 293T cells transfected with a HA-tag protein. Formalin-fixed cells were permeabilized with 0.1% Triton X-100 in TBS for 5-10 minutes and blocked with 3% BSA-PBS for 30 minutes at room temperature. Cells were probed with the primary antibody in 3% BSA-PBS and incubated overnight at 4 °C in a humidified chamber. Cells were washed with PBST and incubated with a DyLight 594-conjugated secondary antibody (red) in PBS at room temperature in the dark. DAPI was used to stain the cell nuclei (blue).