

WB-Validated GNB2 Lentiviral shRNA Knockdown Kit



Catalog #: V63512

Aliases

GNB2; Transducin Beta Chain 2; Signal-Transducing Guanine Nucleotide-Binding Regulatory Protein Beta Subunit 2; Guanine Nucleotide Binding Protein (G Protein), Beta Polypeptide 2; Guanine Nucleotide-Binding Protein G(I)/G(S)/G(T) Beta Subunit 2; Guanine Nucleotide-Binding Protein G(I)/G(S)/G(T) Subunit Beta-2; G Protein, Beta-2 Subunit; G Protein Subunit Beta-2; Heterotrimeric Guanine Nucleotide-Binding Protein 2C1; Epididymis Secretory Sperm Binding Protein; SSS4; NEDHYDF; HG2C1;SSS4

Background

Gene Name: GNB2

NCBI Gene Entry: [2783](#)

Storage

Store at -80 °C for one year.

Kit Components

1. WB-Validated GNB2 shRNA lentiviral particles (4 mL)
2. Non-Target shRNA lentiviral particles (4 mL)
3. Verification Tool: KD-Validated Anti-G protein subunit beta 2 Rabbit mAb #63512 (5 µL)

Tested Cell Line

HeLa

Validation Methods

RT-qPCR; Western Blotting (WB)

Shipping

Shipped with dry ice. Immediately store the product in a standard freezer at -80°C upon receipt.

Instructions For Use

The following protocol uses HeLa cell as an example assuming your cell culture medium is DMEM-based.

1. Release 0.5 million HeLa cells into a 35 mm tissue culture dish in 2 mL of the growth medium (DMEM containing 10% FBS and 1% pen/strep). Cell density should reach 50-60% confluence the following day.
2. 24 h after cell release, pre-warm the shRNA lentiviral medium to 37°C.

SUPPORT

SUPPORT@GENUINBIOTECH.COM
TEL: +1-540-855-7041

ORDERS


SALES@GENUINBIOTECH.COM
FAX: +1-540-855-7041

WWW.GENUINBIOTECH.COM

3. Discard 1 mL of the original growth medium of the 35 mm dish.
4. Using a serological pipette, gently mix the lentiviral solution 3 times.
5. Carefully add 1 mL of the lentiviral solution to the well. Tip: To prevent splashing, add the solution to the dish along the wall.
6. Add a polybrene stock solution to the culture medium at a final concentration of 5 µg/mL. Gently swirl the dish to mix.
7. 48 h after cell release, without discarding the original medium, add another 1 mL of lentiviral medium directly into the dish.
8. Add an additional polybrene stock solution into the dish to obtain a final concentration of 5 µg/mL. Tip: Now, the medium in the dish should be a total of 3 mL.
9. 72 h after cell release, cells may reach confluence. Trypsinize the cells off the 35 mm dish and culture those cells in a 60 mm dish.
10. Add puromycin to the dish at a final concentration of 4 µg/mL. Tip: To assess the efficacy of puromycin selection, culture a dish of wild-type HeLa cells as a negative control.
11. Allow puromycin selection for 48 h. Almost all wild-type HeLa cells should die, while the dish infected with lentiviruses should have some remaining cells.
12. Replace the medium with regular growth medium without puromycin and allow the cells to grow to confluence before harvesting or staining.

- Note:**
1. This product is for research use only.
 2. This product is only supplied to end users.
 3. Do not use this product for commercial.

Validation Data

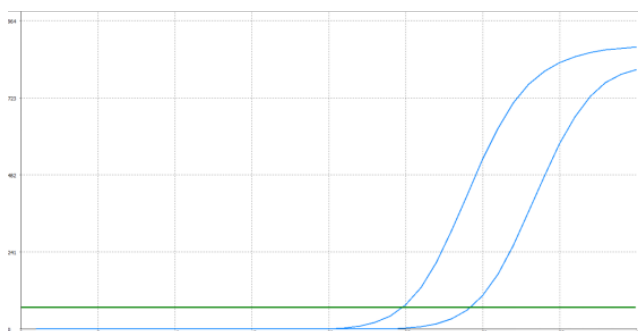

Genotype	Ct Value
Wild-Type	24.63
Knock-Down	28.94
ΔCt (CtKD-CtWT)	4.31
% mRNA Reduction	95%

Copyright ©2025 Genium Biotechnologies LLC

RT-qPCR analysis. HeLa cells were infected with GNB2-specific shRNA lentiviral particles, total RNA was extracted from wild-type and knockdown cells, RT-qPCR was performed using gene-specific primers. ΔC_t (CtKD-CtWT) was used to calculate mRNA reduction (%) between wild-type and knockdown cells using the following formula: $(1 - 1/2^{\Delta C_t}) \times 100\%$.


SUPPORT

SUPPORT@GENUINBIOTECH.COM
TEL: +1-540-855-7041

ORDERS

SALES@GENUINBIOTECH.COM
FAX: +1-540-855-7041

WWW.GENUINBIOTECH.COM

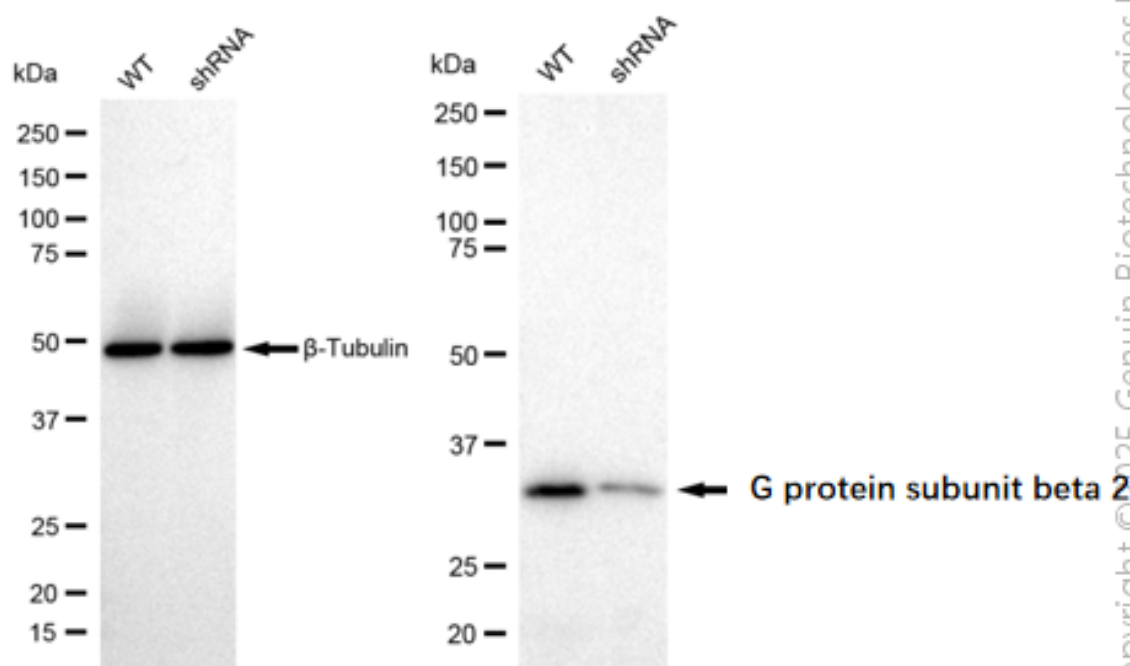

Western blotting analysis. GNB2 protein expression in wild-type (WT) and shRNA knockdown (KD) HeLa cells was detected using Western blotting. β-Tubulin served as a loading control. The blots were incubated with primary antibodies against GNB2 and β-Tubulin, respectively, followed by incubating with HRP-conjugated goat anti-rabbit secondary antibody. Images were developed using FeQ™ ECL Substrate Kit.