

WB-Validated CUL5 Lentiviral shRNA Knockdown Kit



Catalog #: V65564

Aliases

CUL5; Cullin 5; VACM-1; Vasopressin-Activated Calcium-Mobilizing Receptor 1; Cullin-5; CUL-5; VACM1; Cullin-5 (Vasopressin-Activated Calcium-Mobilizing Receptor-1); Vasopressin-Activated Calcium-Mobilizing Receptor-1

Background

Gene Name: CUL5

NCBI Gene Entry: [8065](#)

Storage

Store at -80 °C for one year.

Kit Components

1. WB-Validated CUL5 shRNA lentiviral particles (4 mL)
2. Non-Target shRNA lentiviral particles (4 mL)
3. Verification Tool: KD-Validated Anti-Cullin 5 Rabbit mAb#65564 (5 µL)

Tested Cell Line

293T

Validation Methods

RT-qPCR; Western Blotting (WB)

Shipping

Shipped with dry ice. Immediately store the product in a standard freezer at -80°C upon receipt.

Instructions For Use

The following protocol uses 293T cell as an example assuming your cell culture medium is DMEM-based.

1. Release 0.5 million 293T cells into a 35 mm tissue culture dish in 2 mL of the growth medium (DMEM containing 10% FBS and 1% pen/strep). Cell density should reach 50-60% confluence the following day.
2. 24 h after cell release, pre-warm the shRNA lentiviral medium to 37°C.
3. Discard 1 mL of the original growth medium of the 35 mm dish.
4. Using a serological pipette, gently mix the lentiviral solution 3 times.
5. Carefully add 1 mL of the lentiviral solution to the well. Tip: To prevent splashing, add the

SUPPORT

SUPPORT@GENUINBIOTECH.COM
TEL: +1-540-855-7041

ORDERS

SALES@GENUINBIOTECH.COM
FAX: +1-540-855-7041

WWW.GENUINBIOTECH.COM

solution to the dish along the wall.

6. Add a polybrene stock solution to the culture medium at a final concentration of 5 µg/mL.

Gently swirl the dish to mix.

7. 48 h after cell release, without discarding the original medium, add another 1 mL of lentiviral medium directly into the dish.

8. Add an additional polybrene stock solution into the dish to obtain a final concentration of 5 µg/mL. Tip: Now, the medium in the dish should be a total of 3 mL.

9. 72 h after cell release, cells may reach confluence. Trypsinize the cells off the 35 mm dish and culture those cells in a 60 mm dish.

10. Add puromycin to the dish at a final concentration of 4 µg/mL. Tip: To assess the efficacy of puromycin selection, culture a dish of wild-type 293T cells as a negative control.

11. Allow puromycin selection for 48 h. Almost all wild-type 293T cells should die, while the dish infected with lentiviruses should have some remaining cells.


12. Replace the medium with regular growth medium without puromycin and allow the cells to grow to confluence before harvesting or staining.

Note: 1. This product is for research use only.

2. This product is only supplied to end users.

3. Do not use this product for commercial.

Validation Data

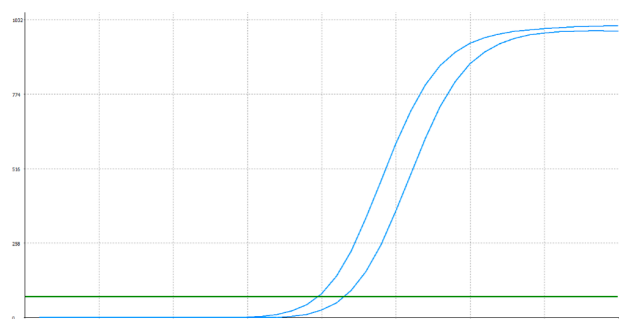


Genotype	Ct Value
Wild-Type	19.72
Knock-Down	21.45
$\Delta C_t (C_{tKD} - C_{tWT})$	1.73
% mRNA Reduction	↓ 69%

Copyright ©2025 Genuin Biotechnologies LLC

RT-qPCR analysis. 293T cells were infected with CUL5-specific shRNA lentiviral particles, total RNA was extracted from wild-type and knockdown cells, RT-qPCR was performed using gene-specific primers. $\Delta C_t (C_{tKD} - C_{tWT})$ was used to calculate mRNA reduction (%) between wild-type and knockdown cells using the following formula: $(1 - 1/2^{\Delta C_t}) \times 100\%$.

WB-Validated CUL5 Lentiviral shRNA Knockdown Kit

PAGE 3

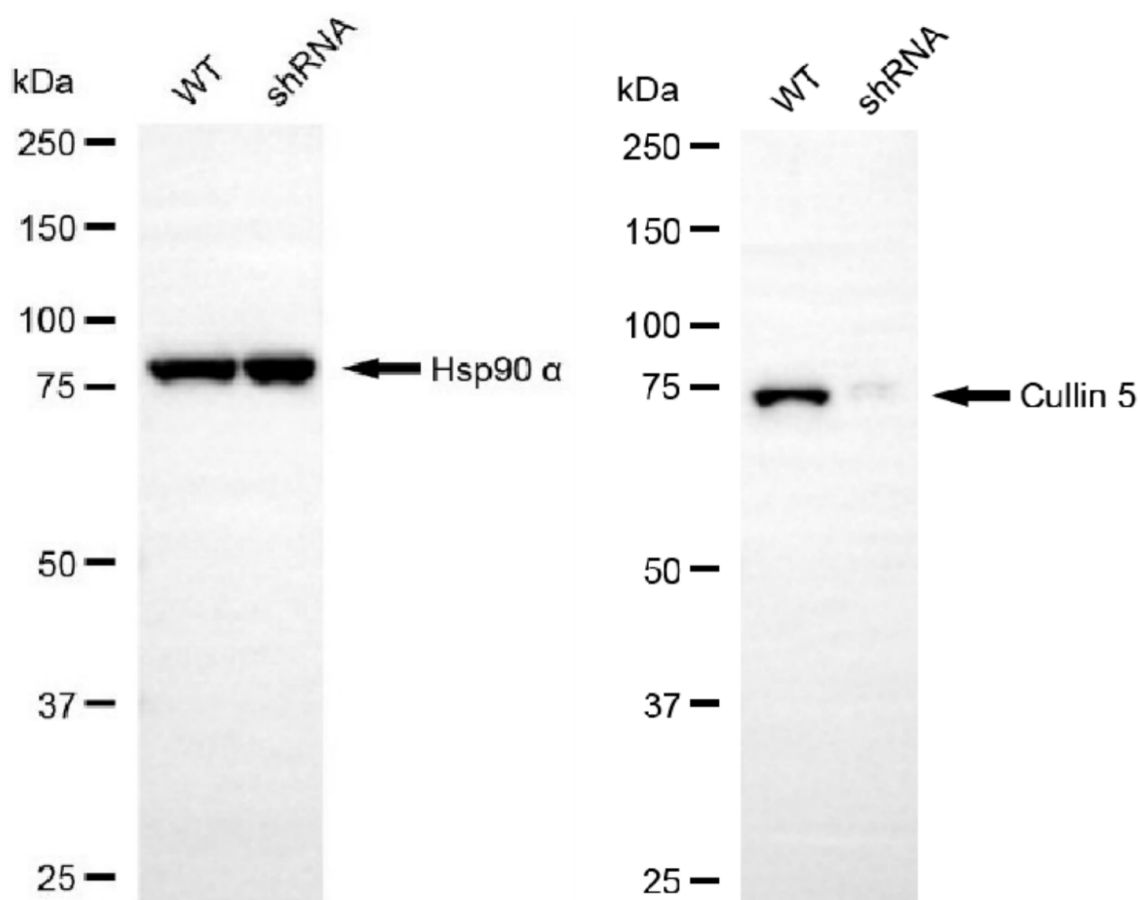

Copyright ©2025 Genuin Biotechnologies LLC

Western blotting analysis. CUL5 protein expression in wild-type (WT) and shRNA knockdown (KD) 293T cells was detected using Western blotting. Hsp90 α served as a loading control. The blots were incubated with primary antibodies against CUL5 and Hsp90 α , respectively, followed by incubating with HRP-conjugated goat anti-rabbit secondary antibody. Images were developed using NaQ™ ECL Substrate Kit.

SUPPORT

SUPPORT@GENUINBIOTECH.COM
TEL: +1-540-855-7041

ORDERS

SALES@GENUINBIOTECH.COM
FAX: +1-540-855-7041

WWW.GENUINBIOTECH.COM